



Now Available with
Masimo SET®



Welch Allyn Propaq® LT

Small enough to fit in the palm of your hand.
Rugged enough to perform in even the
toughest situations.

WelchAllyn®

Advancing Frontline Care™

Welch Allyn Propaq® LT.

The one monitor you can deploy when –
and where – you need it most.

Ease of use, durability and flexibility combine
to give you the perfect monitor for today's
emergency response situations.

Weighing in under 1 kg, the Propaq LT is a
multipurpose patient monitor that's flexible enough
to be used in all types of terrain and environments.
Compact enough to fit in the palm of your hand, the
Propaq LT uses time-tested, durable Propaq technology
to give you the accurate vital signs readings you need.
Plus the Welch Allyn Smartcuf® NIBP technology
provides accurate readings during motion.

And since it's tested to withstand a 1.8 m drop,
it can easily endure all of the bumps and bruises
that are a regular part of your day and still provide
accurate patient vitals information. What's more,
the Propaq LT can be packaged with a Welch Allyn
AED 10™ in a single bag, giving you a powerful
monitor/ defibrillator combination that's easy to carry.

The robust, lightweight Ambulance Mount
gives emergency caregivers fast, easy
access to the device while maintaining
safe, secure transport monitoring.

The lightweight, portable Propaq LT is ideal
for any terrain and environment, even when
you have other equipment to carry.





SHOWN ACTUAL SIZE:
The Welch Allyn Propaq LT has shock resistance of up to 75 G.

Display Modes Alarm/Silence Reset

Full-Disclosure Snapshot Navigational Keypad Blood Pressure

Advanced Masimo features include Signal Identification and Quality Indicator (SIQ) and Perfusion Index (PI) shown below in tabular trend with waveform display.



ID: SET_NEASILJJ Adult Rm 20:56:08

SpO2 2x

Tab/wave	80	20	96
Time	HR/min	NIBP mmHg	Resp/min SpO2 %
20:56	80	n/a	20 96 (5.91)
20:55	80	n/a	20 96 (5.92)
20:54	80	n/a	20 96 (5.92)
20:53	80	n/a	20 96 (5.92)
20:52	80	n/a	20 96 (5.91)
20:51	80	n/a	20 96 (5.90)
20:50	80	n/a	22 96 (5.91)

Quickly toggle between three display modes that are right for your care area (five to choose from). Shown on monitor are large numerics, and above are two waveform and tabular trend with waveform displays.

Use the AED 10 in combination with the Propaq LT to monitor your patients and also be prepared in case of cardiac arrest.



Key Features

- Weighs less than 1 kg
- Shock resistance of up to 75 G – enough to withstand a 1.8 m drop
- Up to 24 hours of battery life
- Heart rate, 3- or 5-lead ECG, SpO₂, NIBP and respiration
- Masimo SET® motion-tolerant SpO₂ or Nellcor® SpO₂
- Smartcuf® motion-tolerant NIBP
- Neonatal, paediatric and adult patient modes
- Continuous monitoring and spot-check modes
- Colour LCD display with backlight, display timeout and lockout controls
- Fast and intuitive camcorder-style controls
- 24-hour tabular and graphical trending with one-minute sampling rate
- Storage of up to 20 full-disclosure data snapshots
- Transfers information to a PC for printing
- Large colour display option

Ambulance Mount

The Starmed ambulance mounting system securely fixes the Propaq LT during transport.

- Lightweight, robust ambulance mounting solution for use in EMS and Military environment
- Provides fast, easy access
- Maintains safe, secure transport monitoring
- Tested and certified for EN 1789 (Shock and Vibration) and CS 27.561 (Emergency Landing Conditions)




A powerful combo

Welch Allyn Propaq LT and Welch Allyn AED 10™ – a powerful monitor/defibrillator combo that's easy to carry in a single bag.



Incorporating the latest technologies

The Propaq LT can communicate with any Windows® 2000/XP PC via a USB port, allowing you to custom-configure the device, continuously monitor a patient, and print tabular trends and snapshots of patient data. Tabular trend and snapshots of patient data can be converted to PDF files for easy integration into electronic medical records (EMRs).

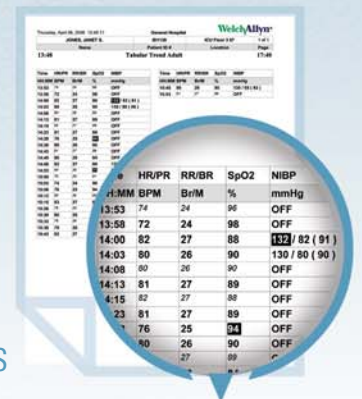

Print an A4 tabular trend report to show patient data, highlight alarm events, and attach to patient records, or print patient snapshot waveforms.

For more information on the Welch Allyn Propaq LT, please contact your local Welch Allyn sales representative or visit us online at www.welchallyn.com.



All Propaq vital signs monitors come with Smartcuf technology, a new standard in NIBP accuracy. This patented software improves NIBP performance and accuracy in the presence of motion artifact.



This monitor is approved for use during all phases of flight aboard U.S. Army aircraft. MIL STD 461E
Army: CE101, CE102, CS101, CS114, CS115, CS116, RE102, RE103



Only Masimo SET SpO₂ sensors, including LNCS series of cables and NR value disposable sensors, should be used with Masimo pulse oximetry. Masimo, Masimo SET, LNCS, and NR are registered trademarks of Masimo Corporation.



Only Nellcor® pulse oximetry sensors, including Oximax® patient-dedicated adhesive sensors, should be used with the Nellcor pulse oximetry option. Nellcor is a registered trademark of Nellcor Puritan Bennett Inc.

UK Customer Services
Welch Allyn UK Ltd, Cublington Road, Aston Abbotts,
Buckinghamshire HP22 4ND, England
tel: +44 (0) 207 365 6780 fax: +44 (0) 207 365 9694

Welch Allyn South Africa (Pty) Ltd, 394 Surrey Avenue,
Ferndale Randburg, 2194, South Africa
P.O. Box 809, Fontainebleau, South Africa 2032
tel: +27 (0) 11 777-7555 fax: +27 (0) 11 777-7556

Welch Allyn European Customer Services
Navan Business Park, Dublin Road, Navan, Co. Meath, Republic of Ireland
tel: +353 46 90 67790 fax: +353 46 90 67754



The Propaq LT Ambulance Mount is manufactured exclusively by Starmed for the Welch Allyn Propaq LT



Advancing Frontline Care™