

Hear more of what you
need to hear and less of
what you don't.

Welch Allyn Master Elite Stethoscope



An innovative sensor-based stethoscope—designed with veterinarians in mind.

Unlike traditional mechanical or electronic stethoscopes, our Master Elite Stethoscope gives you precise sound pickup, eliminating sound loss. With an excellent signal-to-noise ratio and two predefined filter settings, it allows you to pick up those faint crackles, wheezes, soft murmurs, and high-pitched pulmonary sounds that may be difficult to hear with traditional stethoscopes. And the Master Elite Stethoscope provides excellent sound quality even through thick tissue, feathers, and fur.

WelchAllyn[®]

Advancing Frontline Care[™]

Welch Allyn Master Elite Stethoscope


The high-tech stethoscope with excellent sound quality.

With its unmatched clarity of heart and pulmonary sounds, the moisture-resistant Master Elite enables veterinarians to listen, analyze, and communicate sounds with complete confidence. Featuring a two-position frequency selection switch, it allows you to listen using either low (bell) or high (diaphragm) frequencies. In either mode, you can adjust the volume without affecting sound quality.

- Dual-position filter switch lets you listen to high or low frequencies
- Adjustable binaurals and eartips ensure all-day comfort
- Long-life lithium ion battery lasts up to 140 hours
- 22" tube length (4" longer than our human version)
- One-year warranty

SensorTech™—The technology that makes conventional stethoscopes obsolete.

Located in the center of the stethoscope's chestpiece to help reduce ambient noise, the sensor itself is small yet powerful enough to pick up the faintest of signals, even on tiny animals and uneven surfaces. And the stethoscope's predefined frequency filters may help you hear the types of sounds that can enhance your ability to diagnose.

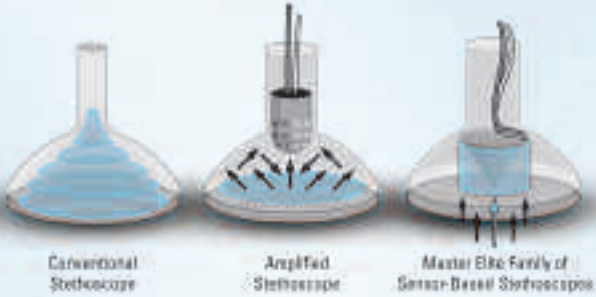

Transducer placement

Master Elite's surface placement eliminates unwanted resonance and vibration.

Welch Allyn Master Elite stethoscopes feature patented sensor technology and an exclusive surface transducer placement that virtually eliminates unwanted resonance and vibration, enhancing your ability to accurately make appropriate diagnoses.

ordering info

Welch Allyn Master Elite Stethoscope (with 22" Tube)	5079-405V
Master Elite Analyzer Software with ECG	5079-430V
Master Elite Distributor	5079-401
Master Elite Headset (for use with the Distributor)	5079-414

specifications

LISTENING MODES	bell/diaphragm	VOLUME CONTROL	15 - 90 dB
BATTERY LIFE	1 year	ON/OFF BUTTON	timer function – 3 minutes
WARRANTY	1 year	DISTRIBUTOR	optional
WEIGHT	170 g	ANALYZER	optional
TUBE LENGTH	22"	CLEANING	cotton/chlorhexidine
		ALLERGY	does not contain latex

For more information call your local representative or Welch Allyn Customer Service at 1.800.535.6663.

4341 State Street Road, PO Box 220, Skaneateles Falls, NY 13153-0220 USA
 (p) 1.800.535.6663 (f) 1.315.685.3361 www.welchallyn.com/dvm