

Hear the Difference.



## Welch Allyn Veterinary Stethoscopes

FOR VETERINARIANS



Hear the difference from the leader in diagnostic care.

For over a decade, we have been focused on providing veterinarians with exceptional performance and value from every Welch Allyn device. Our complete lineup of stethoscopes is no exception.

From the innovative Elite® Electronic sensor-based stethoscope to the Harvey™ DLX with trumpet brass construction, every Welch Allyn product has been designed to provide you with outstanding value and performance.



**WelchAllyn**®


Advancing Frontline Care™



## Elite® Electronic Sensor-based Stethoscope

Unlike traditional amplified stethoscopes that use a microphone in the head to collect the sound, the Elite Electronic places the acoustical sensor directly in contact with the diaphragm, allowing it to pick up even the faintest of sounds.

- 2-position filter switch lets you listen to high or low frequencies
- Adjustable binaurals and ear tips ensure all-day comfort
- Long-life lithium ion battery lasts up to 140 hours
- Fits comfortably around your neck or in your pocket
- Hear through different surfaces (scales)

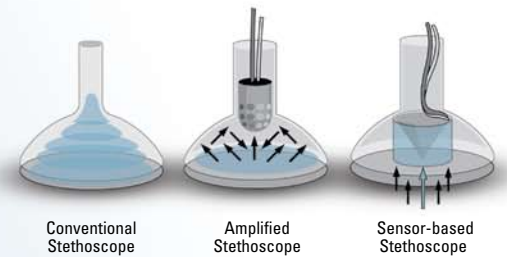

**Transducer Placement**  
Elite Electronic's surface placement eliminates sound loss



## Harvey™ DLX Stethoscope

Trumpet brass construction translates to hearing more of the important cardiac, pulmonary and vascular sounds without relying on difficult "variable pressure" techniques.

- Full-range cardiology/pulmonology stethoscope
- Available in double- or triple-head designs
- Chrome-plated brass chestpiece and binaurals
- Rotatable binaurals and interchangeable comfort-sealing ear tips
- Includes pediatric diaphragm for auscultation of large and small animals



## Harvey™ Elite® Stethoscope

The Harvey Elite offers updated styling with superior comfort and a stainless steel head that transmits sounds crisply and clearly.

- Full-range cardiology/pulmonology stethoscope
- Double head (bell and diaphragm)
- Rotatable binaurals and interchangeable comfort-sealing ear tips

**5079-232 Harvey Deluxe Accessory Kit**  
(Includes one pediatric chestpiece, one flat diaphragm, one corrugated diaphragm, one pair soft eartips, one pair hard eartips)

## Welch Allyn Stethoscope Specifications and Ordering Information

| Product Description | ELECTRONIC STETHOSCOPE         | CARDIOLOGY GRADE   |   |
|---------------------|--------------------------------|--|---|
|                     | ELITE ELECTRONIC               | HARVEY DLX   | HARVEY ELITE  |
| Product Type        | Electronic Stethoscope         | Double/Triple Head   | Double Head   |
| Head Diameter       | 1.5" (38.5 mm)                 | 1.9" (48.3 mm)   | 1.9" (48.3 mm)  |
| Overall Length      | 32.5" (825.5 mm)               | 25" (635 mm)   | 25" (635 mm)  |
| Tubing Length       | 21.5" (546.1 mm)               | 16.5" (419 mm)   | 16.5" (419 mm)  |
| Latex Free          | Yes                            | Yes  | Yes   |
| Chestpiece Material | Plastic                        | Brass  | Stainless Steel   |
| Diaphragm Material  | Plastic                        | Polycarbonate  | Polycarbonate   |
| Diaphragm Style     | Flat                           | Flat or Corrugated   | Flat or Corrugated  |
| Type of Tubing      | Single                         | Dual Lumen   | Dual Lumen  |
| Anti-Chill Ring     | No                             | No   | Yes   |
| Warranty            | 1-year                         | 10-year  | 10-year   |
| Product Numbers     | 5079-405V<br>5079-430 Analyzer | <ul style="list-style-type: none"> <li>● 5079-321 Triple, Black</li> <li>● 5079-322 Triple, Burgundy</li> <li>● 5079-323 Triple, Navy</li> <li>● 5079-324 Triple, Forest Green</li> <li>● 5079-325 Double, Black</li> <li>● 5079-326 Double, Burgundy</li> <li>● 5079-327 Double, Navy</li> <li>● 5079-328 Double, Forest Green</li> </ul> | <ul style="list-style-type: none"> <li>● 5079-125 Elite, Black</li> <li>● 5079-270 Elite, Burgundy</li> <li>● 5079-271 Elite, Navy</li> <li>● 5079-284 Elite, Forest Green</li> </ul> |

To find out how Welch Allyn Stethoscopes can help you hear the difference, go to [www.welchallyn.com/dvm](http://www.welchallyn.com/dvm) or call 1.800.535.6663.

4341 State Street Road, Skaneateles Falls, NY 13153 USA  
(p) 800.535.6663 (f) 315.685.3361 [www.welchallyn.com/dvm](http://www.welchallyn.com/dvm)

**WelchAllyn®**

Advancing Frontline Care™